




Injection Systems



Nipro Needles and Syringes

The challenge of any manufacturer in today's demanding healthcare environment is to provide customers with products that deliver both exceptional performance and value. Nipro lives up to that promise with a full range of injection systems that match a variety of clinical needs. Nipro's state-of-the-art manufacturing assures safe products of consistently high quality.

Disposable Syringe

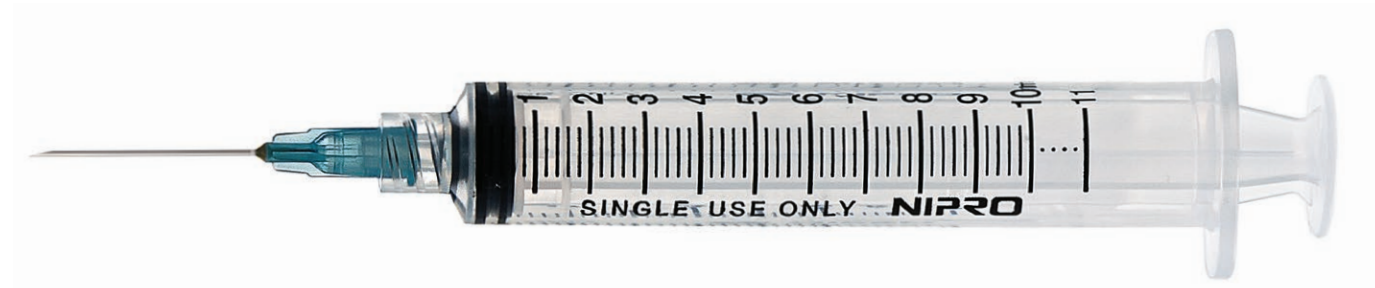

- Polypropylene syringe barrel and specially formulated elastomer stopper result in a smooth, accurate motion along the entire length of the barrel.
- Clear barrel assures rapid, precise volume reading.
- ETO or Gamma sterilized and undergo automated inspection to assure uniform quality.
- Available with Luer Lock, Luer Slip and Eccentric Luer Tip.

SafeTouch™ Safety Needle



- Double-locking mechanism with audible “click” ensures needle is permanently covered.
- Compatible with all manufacturers’ (ISO 594 compliant) syringe hubs.
- Bevel-up needle indication for low-angle injection and blood draws.
- Color-coded needle hub and shield for quick identification of needle gauge.

Hypodermic Needle



- Proprietary grinding process results in an exceptionally sharp lancet point able to penetrate tissue smoothly and with minimal trauma.
- Needle wall is ultra-thin compared to a conventional needle, permitting rapid and efficient drug delivery.
- Translucent hub provides confirmation of flashback.
- Color-coded hub for accurate, quick identification of needle gauge.

Dose-Saver® Syringe



- Low “dead space” at upper end of barrel reduces waste of sometimes costly medications.
- Highly translucent syringe barrel enhances readability.
- Elastomer stopper helps assure consistently smooth plunger action.

ORDERING INFORMATION

Hypodermic Needle

Item #	Description	Packaging
AH+1825	18GX1	
AH+1838	18GX1.1/2	
AH+1925	19GX1	
AH+1938	19GX1.1/2	
AH+2025	20GX1	
AH+2038	20GX1.1/2	
AH+2050	20GX2	
AH+2125	21GX1	
AH+2138	21GX1.1/2	
AH+2219	22GX3/4	
AH+2225	22GX1	
AH+2238	22GX1.1/2	
AH+2325	23GX1	
AH+2338	23GX1.1/2	
AH+2516	25GX5/8	
AH+2525	25GX1	
AH+2538	25GX1.1/2	
AH+2613	26GX1/2	
AH+2713	27GX1/2	
AH+2732	27GX1.1/4	
AH+3013	30GX1/2	

100/Box,
10 Boxes/Case (1,000)

100/Box,
20 Boxes/Case (2,000)

100/Box,
10 Boxes/Case (1,000)

Dose-Saver Syringe

Item #	Description	Packaging
JD+01D2238	1cc 22GX1.1/2	
JD+01D2516	1cc 25GX5/8	
JD+03D2238	3cc 22GX1.1/2	

100/Box,
10 Boxes/Case (1,000)

SafeTouch Safety Needle

Item #	Description	Packaging
SAH+1825	18GX1	
SAH+1838	18GX1-1/2	
SAH+2015	20GX1	
SAH+2038	20GX1-1/2	
SAH+2125	21GX1	
SAH+2138	21GX1-1/2	
SAH+2225	22GX1	
SAH+2238	22GX1-1/2	
SAH+2325	23GX1	
SAH+2516	25GX5/8	
SAH+2525	25GX1	
SAH+2538	25GX1-1/2	
SAH+2713	27GX1/2	
SAH+2716	27GX5/8	
SAH+3013	30GX1/2	

100/Box,
8 Boxes/Case (800)

1cc Syringe

Item #	Description	Packaging
JD+01L	TB Luer Lock without Needle	
JD+01T	TB Luer Slip without Needle	
JD+01T2516	TB 1cc 25GX5/8	
JD+01T2609	TB 1cc 26GX3/8	
JD+01T2713	TB 1cc 27GX1/2	

100/Box,
10 Boxes/Case (1,000)

3cc Syringe

Item #	Description	Packaging
JD+03L	3cc Luer Lock without Needle	
JD+03L2025	3cc Luer Lock 20GX1	
JD+03L2038	3cc Luer Lock 20GX1.1/2	
JD+03L2125	3cc Luer Lock 21GX1	
JD+03L2138	3cc Luer Lock 21GX1.1/2	
JD+03L2219	3cc Luer Lock 22GX3/4	
JD+03L2225	3cc Luer Lock 22GX1	
JD+03L2238	3cc Luer Lock 22GX1.1/2	
JD+03L2325	3cc Luer Lock 23GX1	
JD+03L2332	3cc Luer Lock 23GX1.1/4	
JD+03L2338	3cc Luer Lock 23GX1.1/2	100/Box, 10 Boxes/ Case (1,000)
JD+03L2516	3cc Luer Lock 25GX5/8	
JD+03L2525	3cc Luer Lock 25GX1	
JD+03L2538	3cc Luer Lock 25GX1.1/2	
JD+03S	3cc Luer Slip without Needle	
JD+03S2025	3cc Luer Slip 20GX1	
JD+03S2125	3cc Luer Slip 21GX1	
JD+03S2219	3cc Luer Slip 22GX3/4	
JD+03S2225	3cc Luer Slip 22GX1	
JD+03S2325	3cc Luer Slip 23GX1	
JD+03S2516	3cc Luer Slip 25GX5/8	

5cc Syringe

Item #	Description	Packaging
JD+05L	5cc Luer Lock without Needle	
JD+05L2025	5cc Luer Lock 20GX1	
JD+05L2038	5cc Luer Lock 20GX1.1/2	
JD+05L2125	5cc Luer Lock 21GX1	
JD+05L2138	5cc Luer Lock 21GX1.1/2	100/Box, 10 Boxes/ Case (1,000)
JD+05L2225	5cc Luer Lock 22GX1	
JD+05L2238	5cc Luer Lock 22GX1.1/2	
JD+05L2525	5cc Luer Lock 25GX1	
JD+05S	5cc Luer Slip without Needle	

10cc Syringe

Item #	Description	Packaging
JD+10L2025-WEI	10cc Luer Lock 20GX1	
JD+10L2038-WEI	10cc Luer Lock 20GX1.1/2	
JD+10L2125-WEI	10cc Luer Lock 21GX1	
JD+10L2138-WEI	10cc Luer Lock 21GX1.1/2	100/Box, 9 Boxes/ Case (900)
JD+10L2225-WEI	10cc Luer Lock 22GX1	
JD+10L2238-WEI	10cc Luer Lock 22GX1.1/2	
JD+10L-WEI	10cc Luer Lock without Needle	
JD+10S-WEI	10cc Luer Slip without Needle	

20cc Syringe

Item #	Description	Packaging
JD+20ES	20cc Eccentric Slip without Needle	50/Box, 10 Boxes/ Case (500)
JD+20L	20cc Luer Lock without Needle	

30cc Syringe

Item #	Description	Packaging
JD+30ES-G	30cc Eccentric Slip without Needle - Gamma	50/Box, 8 Boxes/ Case (400)
JD+30L-G	30cc Luer Lock without Needle - Gamma	

50cc Syringe

Item #	Description	Packaging
JD+50ES-G	50cc Eccentric Slip without Needle - Gamma	25/Box, 8 Boxes/ Case (200)
JD+50L-G	50cc Luer Lock without Needle - Gamma	

CAUTION: Federal law restricts these devices to sale by or on the order of a physician.



NIPRO MEDICAL CORPORATION
 200 Crossing Blvd., Bridgewater, NJ 08807
 T: +1 (908) 393-7030, F: +1 (908) 393-7031
www.nipro.com

© 2016 Nipro Medical Corp. All Rights Reserved.
 Dose-Saver is a registered trademark of Nipro Corporation.
 Printed in U.S.A. MS-BRUS-INJ Rev 1